

Easy to setup. Easy to use.

# TV•Ears Voice Clarifying TV Speaker™

Instruction Guide

# TV•EARS®

**tvears.com**  
Facebook.com/tvearsinc  
Twitter: @tvears

**Sales Main Line:** 1-888-883-3277  
**Support:** 1-866-611-9934  
**Fax:** 1-888-958-7899  
**Email:** info@tvears.com

# TV•EARS

## TV•Ears Product Catalog

### TV•Ears Long Range Digital Headphone System

Item Number — 12281 | MSRP \$169<sup>95</sup>

This is our new and improved voice clarifying long range digital headphone system — Works on any TV. *Includes one TV•Ears Headphone, one TV•Ears Long Range Digital Headphone Transmitter and accessories.*



### TV•Ears Foam Ear Tips

Item Number — 40738 | MSRP \$19<sup>95</sup>

One Year supply (4 pairs) of replacement foam ear tips for the TV Ears 5.0 headset. Tips should be replaced every 90 days for optimal performance and hygiene. Ear tips are only available through TV Ears.



### TV•Ears Long Range Headset System

Item Number — 22741 | MSRP \$199<sup>95</sup>

Our new and improved Voice Clarifying Long Range Digital Headset system allows you to enjoy listening to the TV at up to 100ft away from your TV — Works on any TV. *Includes one TV•Ears Long Range Digital Headset, one Cradle Transmitter and accessories.*



### TV•Ears Long Range Battery

Item Number — 55810 | MSRP \$19<sup>95</sup>

Replacement battery for the new improved long range headset. Should last between 4-6 hours depending on volume, charge and use. Only works on the Long Range Headset.



### TV•Ears 5 Year limited warranty:

Welcome to the TV•Ears family and thank you for purchasing our product!

### What are the terms of the warranty:

This 5 year limited warranty covers product failures for a period of 5 years after the original purchase. The warranty is only valid if the product is purchased from TV•Ears directly or from one of the TV•Ears authorized resellers. During the 5 year warranty, If the product fails within the first year of the original purchase date, it will be replaced once at the cost of shipping and handling with the same or newer model. Subsequently, for the remainder of the 5 year warranty period, if the product fails it will be replaced for 50% off MSRP with the same or newer model of equal value. The following items are NOT covered under this warranty: ear tips, batteries, and shipping charges. Any replacement parts purchased at 50% off automatically extends your warranty into a new 5 year limited warranty.

### How to receive warranty service or return product:

To receive warranty service or for returns, call 1-866-611-9934 to speak to our friendly American based customer support representatives. For returns, you will need to receive an RMA number and place it clearly and legibly on the outside of the package otherwise returns are not accepted. All shipping charges are NOT covered by TV•Ears

### To activate your **FREE 5 Year** limited warranty:

**Call:** 1-866-611-9934  
**Go online:** [tvears.com/warranty-registration](https://tvears.com/warranty-registration)  
**Or Mail to:** TV Ears, Inc.  
2701 Via Orange Way, Suite 1  
Spring Valley, CA 91978

Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.  
This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

- (1) this device may not cause harmful interference, and
- (2) this device must accept any interference received, including interference that may cause undesired operation.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses, and can radiate radio frequency energy, and if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.



# Connecting the TV Speaker to the Mini TV Transmitter



MINI TV TRANSMITTER TOP VIEW



REAR VIEW

"POWER" | "OPTICAL" | "COAXIAL" | "ANALOG" PORTS

The TV transmitter is set up with Near Field Technology which will allow it to work seamlessly as soon as both the Transmitter and Speaker are turned on and there's an Audio connection between the transmitter and the TV.

Plug in the power cord to the transmitter and make sure it is close to the TV as it needs to be connected to the TV. Once the transmitter is on, the Blue LED Ring around the TV Ears Logo will illuminate. Then connect the mini TV Transmitter with either an Optical cord (included) or an Analog cord (included) depending on the setup and TV being connected.

In rare cases, a coaxial connection is the only available connection, that cord needs to be purchased separately for the transmitter and TV connection to be established.

To connect the mini TV Transmitter and the speaker, now connect the Speaker at your desired location up to 100 FT from the transmitter.

Once the Speaker is also turned on, it should connect automatically to the Transmitter and you will be able to enjoy the TV on your new Voice Clarifying Speaker™ with its Voice Enhancement Technology™.

*Note: Since the mini TV Transmitter distributes the signal via wireless UHF broadcast, it is best to only use one transmitter per room to avoid interference.*



OPTICAL TO OPTICAL CORD



COAXIAL TO COAXIAL



3.5MM TO 3.5MM CORD

# TV•Ears Voice Clarifying TV Speaker™ Instructions



REAR VIEW

## DC In 12V Charging Port:

Plug in the included power cord to power on the Speaker. The speaker also contains a rechargeable battery, and will be charged when the power cord is connected to the speaker.

## 3.5mm Audio In Connection:

For direct audio connection to a device such as an Ipad or an Iphone, plug in the included 3.5mm cord directly to both the TV Speaker and the audio device being used.

*\*When the power is low, it will be shut down automatically.*



FRONT VIEW

## Volume Dial Control:

### Turning the Speaker On/OFF


- ▶ Turn the dial clockwise to turn the speaker ON. When the speaker is turned ON, the blue LED ring around the volume dial will be illuminated. To turn it off, please turn the dial counterclockwise all the way and the blue LED light turns off.

### Volume Control

- ▶ Turn the dial clockwise to increase the volume
- ▶ Turn the dial counterclockwise to decrease the volume.

# TV•Ears Voice Clarifying TV Speaker™ Additional Controls

## USB 2.0 PORT AND USB-C PORT

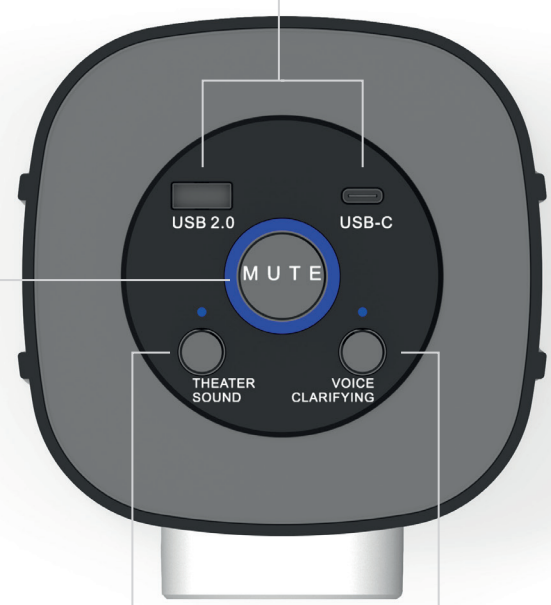

THEATER SOUND

VOICE CLARIFYING

MUTE WITH BLUE LED RING

**USB. 2.0 Port and USB-C Port:** Those two plugs are used to allow charging of other portable devices such as phones or tablets. You will need to use your own cords to charge your devices.

**Mute Button:** To mute the speaker momentarily, press the mute button once. When the blue LED ring is illuminated, that means your speaker is currently muted. Press the mute button again to unmute the speaker.

**Theater Sound Button:** To activate your speaker's theater sound mode for

a better movie experience, press the Theater Sound Button once, you will see the LED light illuminate above the button. To turn the mode off, press the button again.

**Voice Clarifying Button:** To activate your speaker's voice clarifying mode for better dialog recognition, press the Voice Clarifying Button once, you will see the LED light illuminate above the button. To turn the mode off, press the button again.